aquapeel & booster

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Contents

1. To start with	16
1) Overview of product	
2) Product Composition	
3) Detailed Name	
2. How to operate	17
1) Users and areas Not suitable for using product	
2) Button Display	
3) How to install a waste water bin	
4) How to change a Tip	
5) Description of the Tips	
6) How to use each mode	
7) Charging Method	
3. Caution	24
4. Product Specification	25
5. How to solve problems	26
6. How to dispose the product	27
7. Warranty Certificate	28

To start with


Overview of product

Aqua Peel & Booster is a skin care device that has been developed to help users enjoy treatments and massages anytime anywhere that are available only in skin care shops. The peeling mode removes impurities and sebum that were difficult to remove with cleansings to keep your skin clean, and the moisturizing mode provides moisture and nutrition to your skin. If you're too busy to visit skin care shops, Aqua Peel & Booster is the answer to keep your skin clean and radiant, since the device is easy and simple to use.

Product Composition

AQUAPEEL & BOOSTER Product body / Cosmetic solutionX2 (Peeling solutionX1, Moisturizing solutionX1) / Product cleaning solutionX2 / Skin care tipX3 (Face peeling tipX1, T-zone peeling tipX1, Moisturizing tipX1) / Product Cleaning tipX1 / Waste water bin / User manual / Charging adapter

Detailed Name

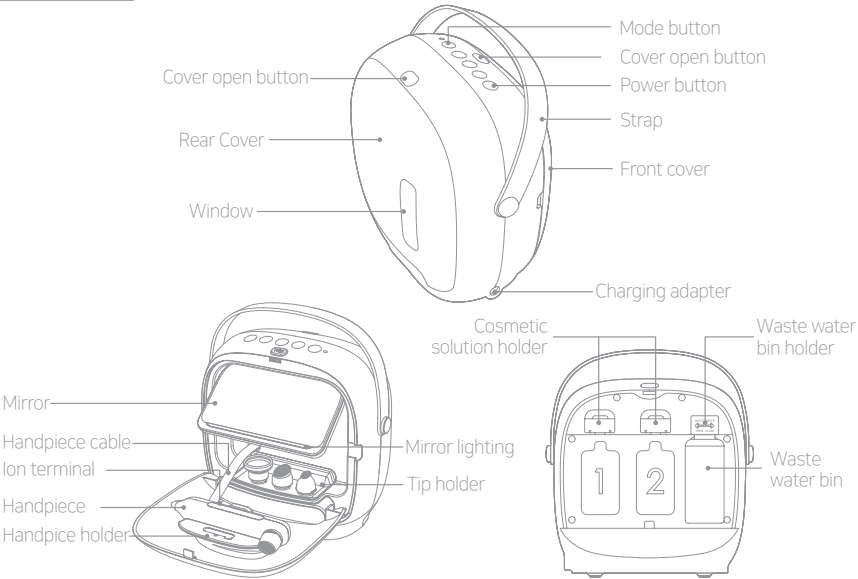How to operate

Users and areas Not suitable for using product

- Excessive tanned skin
- Oral(Uvula) and genital area
- Bleeding area
- Wounded area
- Dermatitis area
- Infected pimples
- Metal allergy or pimples Dermatitis disease
- Patient with major heart diseases.
- Pregnant woman
- Please consult with a doctor in case that user had the skin or plastic surgery. (Including botox and filler procedures)

Button Display

1. Power ON/OFF

Pressing the power button turns on the product (The startup sound is given, the white light illuminates, and the power is turned on). Pressing the power button during operation or in the standby mode turns off the product (The stop sound is given, and the power is turned off).

2. Mode Operation

When powered on, press the required mode button to operate. When the Waste water bin is full, the Water full indicator lights up in red and the operation of Mode stops.

3. Mirror light

You can turn on and off the white LED on the bottom of the mirror with the light button.

4. Water Full Indication Lamp


If the Waste Water Bin is full while the product is in use, the product stops working with an audible signal and red light on the Water full indication lamp. In case that, please operate the mode again after draining the inside of the Waste water bin

How to install a waste water bin

- 1) Press the push button at the top of the back of the product to open the rear cover.
- 2) With the solution handle holding push the container in until it is completely fixed at the holder at the back of the device. Install peeling solution at the 1st holder and the moisturizing solution at the 2nd holder. (Make sure the number indicated on the solution and the number on the back matches)
- 3) Place the Waste water bin on its holder and rotate it to the right until it is completely seated.

* The product does not work if the Waste water bin is not installed correctly.

* If you use other cosmetics other than the solution, it may cause the device to break.

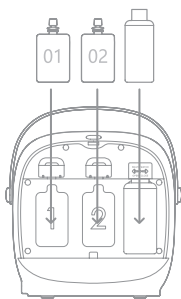Solution

When inserting the solution into the stand, allow the solution handle to hold and push it until it is fully seated.


Waste water bin

Insert the Waste water bin into the groove of the holder and rotate it to the right until you hear a "Click" sound.

How to change a Tip

Replace the tip by pushing and pulling the tip as shown in the image. Use the proper tip for each mode.

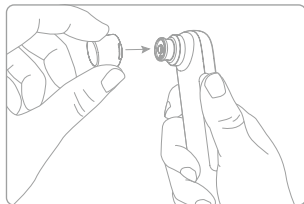Face Peeling Tip

T-zone Peeling Tip

Moisturizing Tip


Product Cleaning Tip

Description of the Tips


1. Face Peeling Tip

In the Peeling mode, the suction tip is used on the face to effectively remove residues on the skin, such as black heads, dead skin cells, and fine dust, to keep the skin clean and smooth.

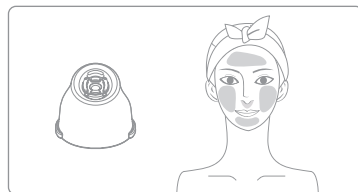
2. T-zone Peeling Tip

In the peeling mode, the T-zone Tip is used to remove blackheads and dead skin cells around the nose effectively.

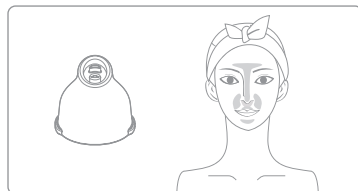
3. Moisturizing Tip

In the Moisturizing mode, the Face Tip is used to help the active ingredients of the cosmetics penetrate deep into the skin through iontophoresis.

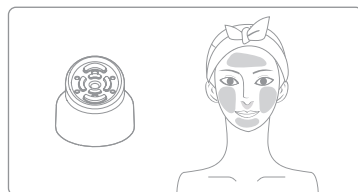4. Product Cleaning Tip

In the Cleaning mode, the Tip is used to clean inside the hose of the product to prevent clogging.

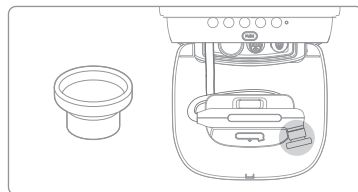
How to use each mode

1. Peeling mode

In this mode, the designated solution and pressure inhaler completely washes away impurities inside your skin that usual cleansing could not remove away.

- 1) Check whether the solution and waste water container are installed on the back of the device properly.
- 2) Insert Face (or T-zone) Peeling Tip into the the hand piece.
- 3) Select the suction level by pressing the Peeling mode button (In a total of three levels) depending on your skin condition (Level 1 for low suction force; Level 2 for medium suction force; and Level 3 for high suction force)
- 4) Choose the pressure level according to your skin condition. There are total 3 pressure levels and the level goes up every time you press the mode button. The pressure level is indicated with a beep sound. (ex: if the beep rings twice when you pressed the button, it means you're on level 2.)
- 5) Hold the handpiece and touch the tip to the back of your hand, palm, or flat, clean skin area, then wait a moment for the solution to come out through the hose.
- 6) When solution starts to come out, massage your skin by attaching and detaching the hand piece on your skin. (If you attach the hand piece on the same area for more than 3 seconds, you may harm your skin due to excessive pressure.)
- 7) We recommend you to use face peeling tip for large areas such as cheek, forehead, chin, and use T-zone peeling tip for curvy areas such as nose.
- 8) The time is set for 10 minutes. When 5 minutes passes the beep rings three times to inform you that half of the total time has passed. (If your skin is sensitive, shorten the time.)

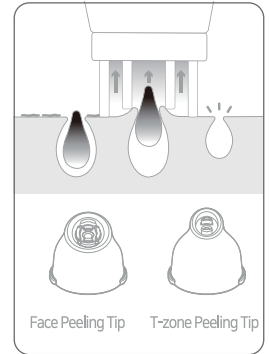
* You can remove waste more effectively by using a face steamer or warm towel to open the pores before using the product.

* Be careful not to use the product in the same area for too long as it may cause irritation or bruising.

2. Moisturizing mode

Micro-current activates your inner skin and provides the solution that contains moisture and nutrition through galvanic feature.

- 1) Insert the Moistuizing Tip into the handpiece.
- 2) While the power is turned on, press the Moistuizing Mode button.
- 3) Press moisturizing mode button while the power is on.
- 4) Choose the pressure level according to your skin condition. There are total 3 pressure levels and the level goes up every time you press the mode button. The pressure level is indicated with a beep sound. (ex: if the beep rings twice when you pressed the button, it means you're on level 2.)
- 5) Hold the handpiece and touch the tip to the back of your hand, palm, or flat, clean skin area, then wait a moment for the solution to come out through the hose.
- 6) When solution starts to come out, massage your skin by attaching and detaching the hand piece on your skin. (if you attach the hand piece on the same area for more than 3 seconds, you may harm your skin due to excessive pressure.)
- 7) The time is set for 5 minutes. After using the device, gently massage your skin to absorb the remaining solution. (If your skin is sensitive, shorten the time.)
- 8) The time is set for 5 minutes. When 3 minutes passes the beep sound rings three times to inform you that half of the total time has passed. (If your skin is sensitive, shorten the time.)

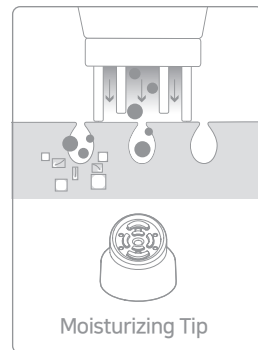
* If you attach the hand piece on the same area for more than 3 seconds, you may harm your skin.


3. Cleaning mode

Clean the product with a cleaning liquid to prevent clogging or bacterial growth inside the product, which may result from residual cosmetics.

- 1) Insert the 2 cleaning solution into holder 1 and 2 at the back of the device and install the Waste water bin.
- 2) Insert the device cleaning tip into hand piece.
- 3) Press cleaning mode button while the power is on.
- 4) The cleaning starts with the beep sound and stops automatically after 3 minutes.
- 5) After cleaning, discard the cleaning solution inside the Waste water bin and wash the container with running water.
- 6) Clean the tips and inside the waste water container with a cotton pad or cotton swab as shown in the images **6**.

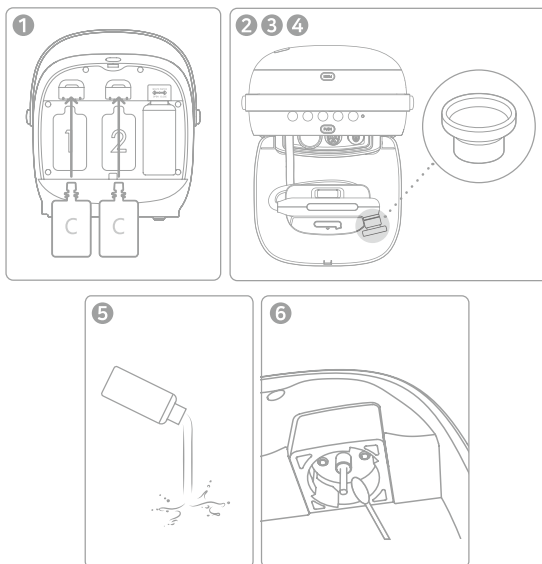
Charging Method

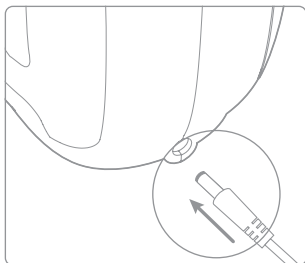
Connect the charging adapter

Insert the charging adapter into the charging connector at the bottom of the side of the product

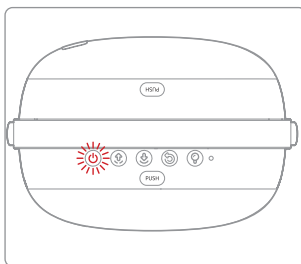Charging in progress

If the charging adapter is inserted correctly, a red light illuminates in the power button.
(The product is available for use during charging.)

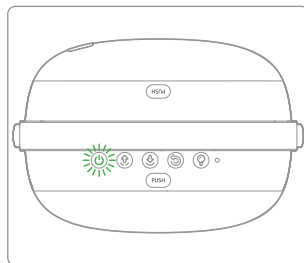Charging complete

When the product is fully charged, a green light illuminates in the power button.

TIP

- * It is recommended to charge using the adapter included in the box. (Charging Power : D.C. 9.0V / 2.0A)
- * If the battery is low during use, a red light illuminates in the power supply and an audible signal sounds. When the battery is completely discharged, the product automatically stops.

Caution

- Understand how the Product operates and use it in the correct sequence.
If you detect any problems on your skin after using the device, stop using the device immediately
- Do not use the Product for any other purpose than its intended use.
- It is recommended to use the device 1-2 times every week. (every time within 20 minutes) the recommended usage time may differ depending on your skin type or condition. Increase or decrease the usage time depending on your skin condition.
- In order to avoid irritation, make sure not to rub or attach on the same area for too long
- Do not plug or unplug with wet hands since it may cause electric shock.
- If you feel any problem during usage, stop using the device immediately, turn off the power, and call the A/S center.
- Do not clean the device with chemicals.
- When not in use, make sure to turn off the power.
- If the product is not used for a long period of time, the battery may completely discharge and the life span of the battery may shorten.
- Do not keep the product on places with high temperature, dusty places, near the fire, or under direct sunlight, since it may cause breakage or shorten the life span of the product.
- If the product breaks, stop using it immediately. If you continue using it, it may cause secondary breakage or cause harm on your skin.
- Make sure the hand pieces do not get tangled during usage. (if you use modes when the cable is tangled, the output voltage may decrease)
- If the cord is damaged, it must be replaced by manufacturer or service department or any similar departments in order to avoid danger.
- When you register the device for A/S, please send it to us after discarding waste

Product Specification

Model Name	COZCORE AQUA PEEL & BOOSTER 2 (MR-CL-213A)
Manufacturer	MIRANG CO.,LTD.
Address	#710, Ace Techno Tower 5 Cha, 20, Digital-ro 31-gil, Guro-gu, Seoul, 08380, Korea
Country of Origin	Republic of Korea
Dimensions (WxDxH)	214x143x225(mm)
Weight	1,421g
Material	A.B.S, PC
Power Consumption	8.5W
Charging Power	DC 9.0V / 2A
Power Supply	Li-ion Battery / 2,000mAh
Manufacturing Date	Marked on the product separately
Customer Service	1833-8393
Registration Number	R-R-MIY-MR-CL-213A
Instructions	As indicated in the user manual

* Battery specification can be changed to other similar one without prior notification.

How to solve problems

Q The pressure is weak while using the device.

A When inserting the solution into the stand, allow the solution handle to hold and press it until it is fully seated.

Q The product does not work.

A It is probable that the battery is discharged. First, charge the battery, then try it again, and if the problem persists, contact the customer center.

Q The product is powered but does not work.

A Make sure the waste water bin is installed correctly. If the waste water bin is not installed correctly, the pressure drops or the product does not operate normally.

Q The product turns off automatically during use.

A The battery is low. Charge the battery and try again.

Q The color in the charging mode

A A red light illuminates during charging and a green light illuminates when charging is complete.

* When the battery is discharged, a red light flickers.

Q The light in the charging mode is not in red or green, also, the product malfunctions.

A Use the adapter provided in the package and check that the product is charged correctly. Other adapters may not charge properly or damage the device. If the problem persists even with the adapter provided, contact the customer center.

Q My skin itches after using the product.

A The product generates micro-current during operation and may cause itching on the skin due to the activation of micro-current.

Q Can I use the Product after a plastic surgery (Botox injection, filler infection, etc.)?

A Consult your physician before use.

Q The area of use has turned red. Is this a problem?


A It is caused by the Product's suction force, and if the problem persists for a long time, please stop using it and consult with doctor.

How to dispose the product

- This product is not household waste.
- When disposing the product, separate batteries from the product and recycle them. It will save natural resources, public health and the environment.
- If you are not sure about how to dispose the product, please contact your local administrative center (resident center).

How to remove the battery (if you're not an expert do not disassemble it)

- 1) Remove the battery only when the power is off and there is no current.
- 2) Separate the front case and the rear case by removing 6 screws on the case of rear side using a driver.
- 3) Remove the 8 connection terminals on the PCB of upper part of the front case, and you will be able to see the battery and the battery holder in the bottom of the body.
- 4) Remove the battery after removing 8 screws on the battery holder with a driver.
- 5) Dispose the removed batteries at a nearby waste battery collection facility.

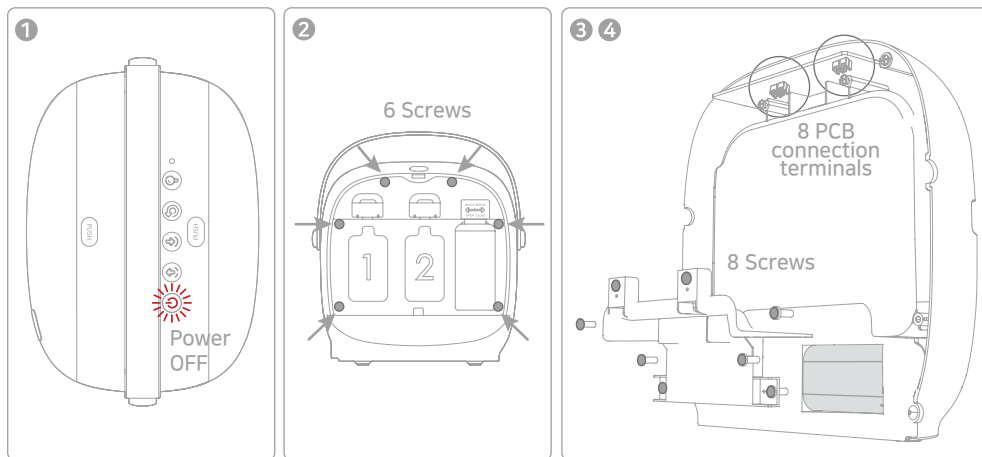Warranty Certificate

- 1. The free warranty of this product will be available only for the failure or defect of the product within 12 months period from the date of purchase. (Battery is 6 month warranty)
- 2. You are required to fill out the date of purchase in this warranty certificate as free warranty period shall be counted based on the date. (In case the date of purchase is not filled out, the corresponding repair cost will be charged even in warranty period)
- 3. Please keep this warranty certificate together with the user manual in a safe place as this certificate will not be re-issued.
- 4. For the following cases, the warranty does not cover the service charge for the repair of the product.
 - 1) Failure caused by negligence of the user in complying with the instructions and caution of the manual in use
 - 2) Failure caused by unauthorized disassembly of this product and external factors such as shock, etc
 - 3) Failure caused by a natural disaster
 - 4) In case the life of consumable goods wore out
 - 5) In case the warranty certificate is not provided
- 5. In case of supply under a separate contract, the warranty is applied according to the contract.

Model Name	COZCORE AQUA PEEL & BOOSTER 2 (MR-CL-213A)
Serial Number	
Customer Name	
Customer Address	
Customer's telephone number	
Date of Purchase	
Purchased Place	

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[Manufacturer] MIRANG Co., Ltd